




OSX Optical Switch

User Manual

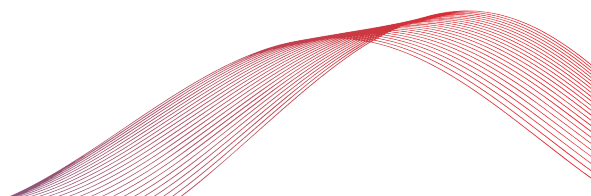

All information contained herein is believed to be accurate and is subject to change without notice. No responsibility is assumed for its use. © JGR Optics Inc., 2020. All rights reserved.

TABLE OF CONTENTS

LIST OF FIGURES AND TABLES.....	v
COMPLIANCE.....	1
FDA-CDRH Compliance.....	1
CSA / IEC Compliance.....	1
CE Compliance.....	1
GENERAL INFORMATION.....	2
OSX Optical Switch Overview.....	2
Applications.....	2
Key Features.....	2
Included Accessories.....	3
SAFETY INFORMATION.....	4
Safety Markings on the Unit.....	4
Classification.....	4
Important Safety Information.....	5
<i>Laser Hazards</i>	5
<i>Electrical Hazards</i>	5
GETTING STARTED.....	7
Initial Inspection.....	7
Operational Requirements.....	7
Product Overview.....	8
<i>OSX Front Panel</i>	8
<i>OSX Rear Panel</i>	8
OPERATION.....	10
Powering Up the Switch.....	10
Using with an RLM.....	10
Establishing Communication.....	10
<i>USB</i>	10
<i>Ethernet</i>	10
OSX WEBPAGE.....	12
Dashboard.....	12
Settings.....	13
<i>Network Settings</i>	13
<i>Upgrade</i>	13
Help.....	14
About.....	15
PROGRAMMING GUIDE.....	17
General Information.....	17
Step-by-step Guide.....	17
Note.....	18

Commands Lists.....	18
MAINTENANCE	21
Cleaning the Unit.....	21
Cleaning the Output Connectors.....	21
STORAGE AND SHIPPING	22
Returning Instruments to JGR Optics	22
Contact Information	22
SPECIFICATIONS	23

LIST OF FIGURES AND TABLES

Figure 1: OSX Optical Switch.....	2
Figure 2: Front view of a 24 channel OSX	8
Figure 3: Rear view of an OSX	9
Figure 4: Retrieving the OSX's IP address	11
Figure 5: OSX webpage – <i>Dashboard</i> tab	12
Figure 6: OSX webpage – <i>Network Settings</i> tab.....	13
Figure 7: OSX webpage – <i>Upgrade</i> tab	14
Figure 8: OSX webpage – <i>Help</i> tab	15
Figure 9: OSX webpage – <i>About</i> tab	16
Table 1: Safety symbols.....	4
Table 2: Environmental requirements	8
Table 3: Detailed description of the OSX rear panel components (see Figure 3).....	9
Table 4: Standard SCPI required commands list.....	19
Table 5: OSX commands list.....	19
Table 6: OSX mechanical and environmental specifications sheet.....	23

1

COMPLIANCE

FDA-CDRH Compliance

Under the US Food and Drug Administration (FDA) Center for Devices and Radiological Health (CDRH), the unit complies with the Code of Federal Regulations (CFR), Title 21, Subchapter J, which pertains to laser safety and labeling. See following link for more information:

- <http://www.accessdata.fda.gov/scripts/cdrh/cfdocs/cfcfr/CFRSearch.cfm?CFRPartFrom=1000&CFRPartTo=1050>

CSA / IEC Compliance

The unit complies with certain standards of the Canadian Standards Association (CSA) and the International Electrotechnical Commission (IEC).

The unit falls in the Installation Category (Overvoltage Category) II under IEC 664. IEC 664 relates to impulse voltage levels and insulation coordination. The category is defined as: local level, appliances, portable equipment, etc., with smaller transient overvoltages than Installation Category (Overvoltage Category) III.

The unit falls in the Pollution Degree 2 category under IEC 1010-1 and CAN/CSA-C22.2 No. 1010.1. The IEC standard on Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use relates to insulation coordination. The CSA standard is on Safety Requirements for Electrical Equipment for Measurement Control, and Laboratory Use, Part I: General Requirements. The Pollution Degree 2 category is defined as follows: “Normally only non-conductive pollution occurs. Occasionally, however, a temporary conductivity caused by condensation must be expected.”

CE Compliance

Electronic test equipment is subject to the EMC Directive in the European Union. The EN61326 standard prescribes both emission and immunity requirements for laboratory, measurement, and control equipment. This unit has undergone extensive testing according to the European Union Directive and Standards.

2

GENERAL INFORMATION

OSX Optical Switch Overview

The OSX Optical Switch (Figure 1) is a customizable benchtop/rackmount instrument, ideal for high-volume manufacturing production testing.

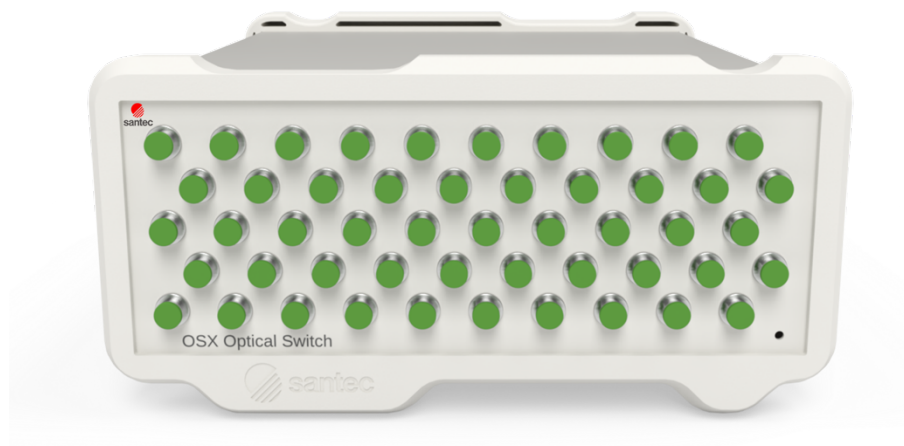


Figure 1: OSX Optical Switch

Available in single-mode or multimode, the switch configurations range from 1x2 up to 1x360. Optional 2x1 or 4x1 built-in switches can provide additional common inputs. The switch may be controlled through remote interface (USB or Ethernet).

Applications

- High volume production testing
- Multiple device testing
- Test automation
- Paired with the RLM for multi-fiber testing

Key Features

- Low IL < 0.7 dB
- ± 0.005 dB repeatability
- Up to 360 outputs

- Chassis modular
- Customizable

Included Accessories

- Test report
- USB A to USB B cable (1.8m)
- Ethernet cable (1.8m)
- AC power cord (2m)

3

SAFETY INFORMATION


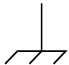

To avoid situations that could result in serious injuries or death, always observe the following precautions.

The safety instructions must be observed whenever the unit is operated, serviced, or repaired. Failure to comply with any of these instructions or with any precaution or warning contained in the user manual is in direct violation of the standards of design, manufacturing, and intended use of the unit. JGR Optics Inc. assumes no liability for the customer's failure to comply with any of these safety requirements.

Safety Markings on the Unit

See Table 1 for symbols and messages that can be marked on the unit. Observe all safety instructions that are associated with a symbol.

Table 1: Safety symbols

	Laser radiation may be present. Refer to the user manual for instructions on handling and operating the unit safely. Avoid looking into any ports near which this symbol appears.
	Frame or chassis terminal for electrical grounding within the unit.
	Protective conductor terminal for electrical grounding to the earth.
WARNING	Procedure can result in serious injury or loss of life if not carried out in proper compliance with all safety instructions. Ensure that all conditions necessary for safe handling and operation are met before proceeding.
CAUTION	Procedure can result in serious damage to or destruction of the unit if not carried out in compliance with all instructions for proper use. Ensure that all conditions necessary for safe handling and operation are met before proceeding.

Classification

The OSX consists of an exposed metal chassis that is connected directly to earth via a power cord and is therefore classified as a Class 1 instrument.

Important Safety Information

Laser Hazards

Warning



- Never look directly into the end of an optical cable connected to an optical output device that is operating. Laser radiation is invisible and direct exposure can severely injure the human eye.

Electrical Hazards

Warning



- Some of the circuits are powered whenever the unit is connected to the AC power source (line power). To ensure that all circuits are powered off, disconnect the power cord from either the power inlet on the unit's rear panel or from the AC line-power source (receptacle). The power cord must always be accessible from one of these points. If the unit is installed in a cabinet, the operator must be able to disconnect the unit from the line power by the system's line-power switch.
- Use only the type of power cord supplied with the unit. If you need to replace a lost or damaged cord, make sure to replace with a power cord of the same type.
- Connect the power cord only to a power outlet equipped with a protective earth contact. Never connect to an extension cord or any receptacle that is not equipped with this feature.
- If using a voltage-reducing autotransformer to power the unit, ensure that the common terminal connects to the earthed pole of the power source.
- Do not interrupt the protective earth grounding. Such action can lead to a potential shock hazard that can result in serious personal injury. Do not operate the unit if an interruption to the protective grounding is suspected.
- Do not operate the unit when its cover or panels have been removed.
- To prevent potential fire or shock hazard, do not expose the unit to any source of excessive moisture.
- Do not use the unit outdoor.
- Operating the unit in the presence of flammable gases or fumes is extremely hazardous.

- If the equipment is used in a manner not specified by the manufacturer, the protection provided by the equipment may be impaired. Only technicians authorized by JGR Optics Inc. should carry out repairs. In addition to voiding the warranty, opening the unit (even when unplugged) can expose you to potential shock hazards.
- Some of the unit's capacitors can be charged even when the unit is not connected to the power source.
- Do not perform any operating or maintenance procedure that is not described in the user manual.

4

GETTING STARTED

Caution



- To avoid injury or death, always observe the precautions listed in SAFETY INFORMATION on page 4.

This manual contains complete operating instructions for safe and effective operation of the OSX Optical Switch. It is recommended that users of the OSX familiarize themselves with contents of this manual before using the instrument.

The inspection report and a description of any customer-requested information may be found in the calibration document envelope included with the instrument.

Initial Inspection

Warning



- To avoid electrical shock, do not initialize or operate the unit if it bears any sign of damage. Ensure that the unit and any devices or cords connected to it are properly grounded.

- ✓ Inspect the package and contents for signs of damage.
- ✓ Ensure all contents are included.
- ✓ Read the user manual thoroughly and become familiar with all safety symbols and instructions to ensure that the unit is operated and maintained safely.
- ✓ If the initial inspection reveals any damage or missing components, immediately notify JGR Optics Inc. and if necessary, the carrier.

Operational Requirements

For the unit to meet the warranted specifications, the operating environment must meet the conditions outlined in Table 2.

Table 2: Environmental requirements

Parameter	Specification
Altitude	Up to 2000m
Temperature	0 to 40°C
Humidity	Up to 95% RH (0 to 40°C)
Voltage	Main supply voltage fluctuations must not exceed $\pm 10\%$ of the nominal voltage

Product Overview

OSX Front Panel

A front view of the OSX switch is shown in Figure 2. The standard front panel bulkhead types are FC, SC or LC.

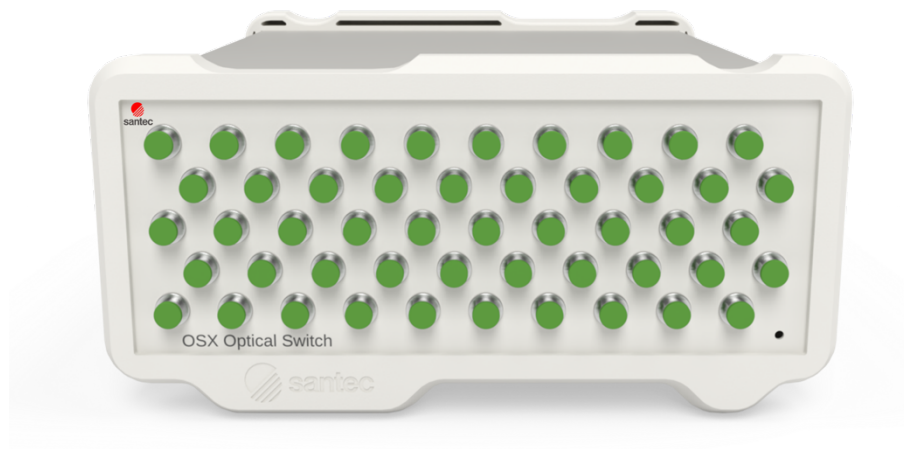


Figure 2: Front view of a 24 channel OSX

OSX Rear Panel

A rear view of the OSX switch is shown in Figure 3. See Table 3 for a detailed description.



Figure 3: Rear view of an OSX

Table 3: Detailed description of the OSX rear panel components (see Figure 3)

Item #	Description
1	LAN/Reset <ul style="list-style-type: none"> Press once: reset network settings
2	Ethernet port <ul style="list-style-type: none"> Connection to LAN
3	USB B port <ul style="list-style-type: none"> Connection to PC or RLM
4	USB A ports <ul style="list-style-type: none"> Connections to peripherals
5	IO switch <ul style="list-style-type: none"> On/off toggle
6	Power input <ul style="list-style-type: none"> Contains user-replaceable fuse

5

OPERATION

Powering Up the Switch

To power up the switch:

1. Verify that the power switch is set to the “off” position (O). Connect the switch to an AC power supply using the power cord provided.
2. Toggle the power switch to the “on” position (I).

Using with an RLM

When using the OSX with an RLM, connect a USB cable from the OSX USB B port to one of the RLM's USB A ports. All switch operations are then controlled via the RLM.

Establishing Communication

If not using with an RLM, the OSX's control is done via USB or Ethernet.

USB

Connect a USB cable from the OSX USB B port to a PC or tablet. VISA drivers are required.

- Recommended: *Rohde and Schwarz*
 - https://www.rohde-schwarz.com/gr/applications/r-s-visa-application-note_56280-148812.html
- Alternatives: *National Instruments, Keysight, etc.*

Ethernet

Each OSX is factory pre-set to use DHCP. To retrieve the IP address of the OSX, it must initially be connected to a PC or tablet via USB and the OSX queried via remote access (see Figure 4):

- Connect the OSX to the network via an Ethernet cable
- Connect the OSX to a PC or tablet with VISA drivers (see above)
- Download, install and run XN1 Terminal:
 - <https://jgropics.com/downloads/software>
- Click on *Detect Devices* to establish communications
- To confirm communications, send the following command:
 - *idn?
- To retrieve the IP address, send the following command:
 - syst:comm:lan:addr?

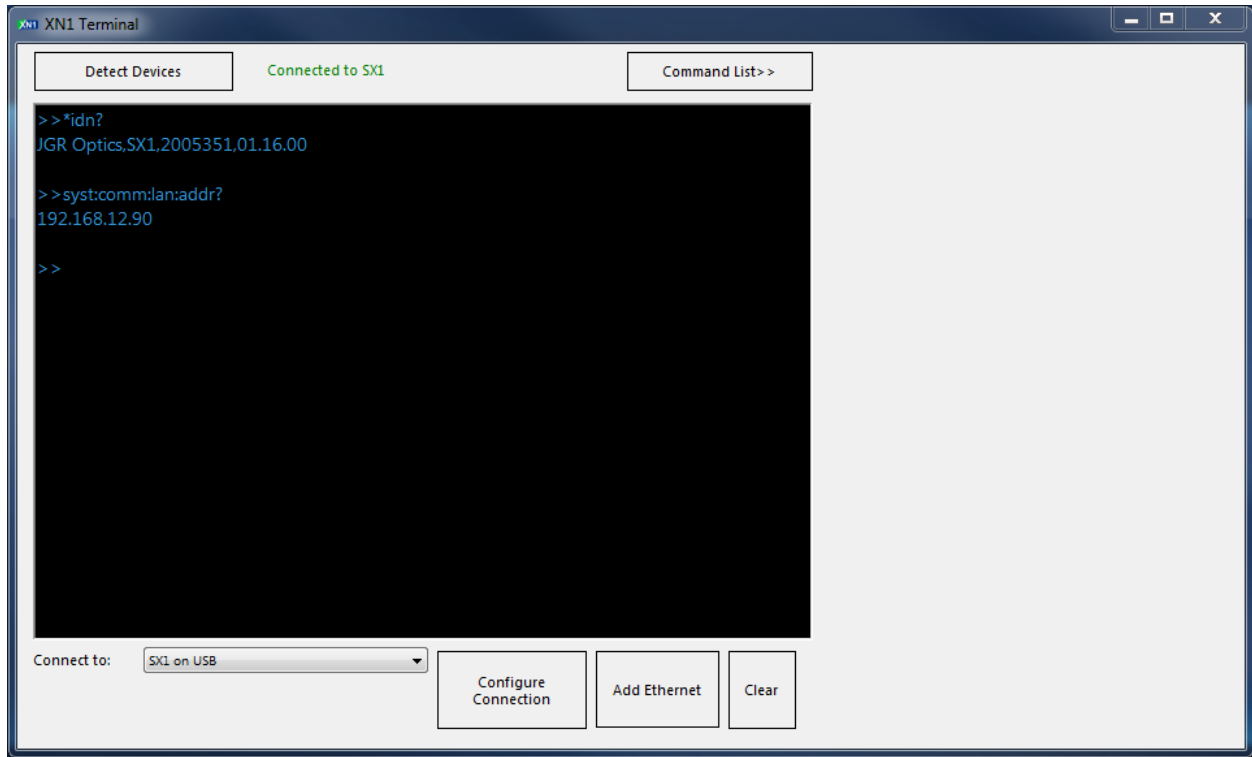


Figure 4: Retrieving the OSX's IP address

6

OSX Webpage

To access the OSX webpage, connect the switch to a network and on any computer or tablet on the same network, open a web browser (recommended: *Google Chrome* or *Firefox*) and enter the OSX's IP address in the URL bar. For more information on retrieving the OSX's IP address, see page 10.

Dashboard

Access the *Dashboard* tab to view the manufactured date, part and serial numbers (Figure 5). You can also download a PDF of the switch test report.

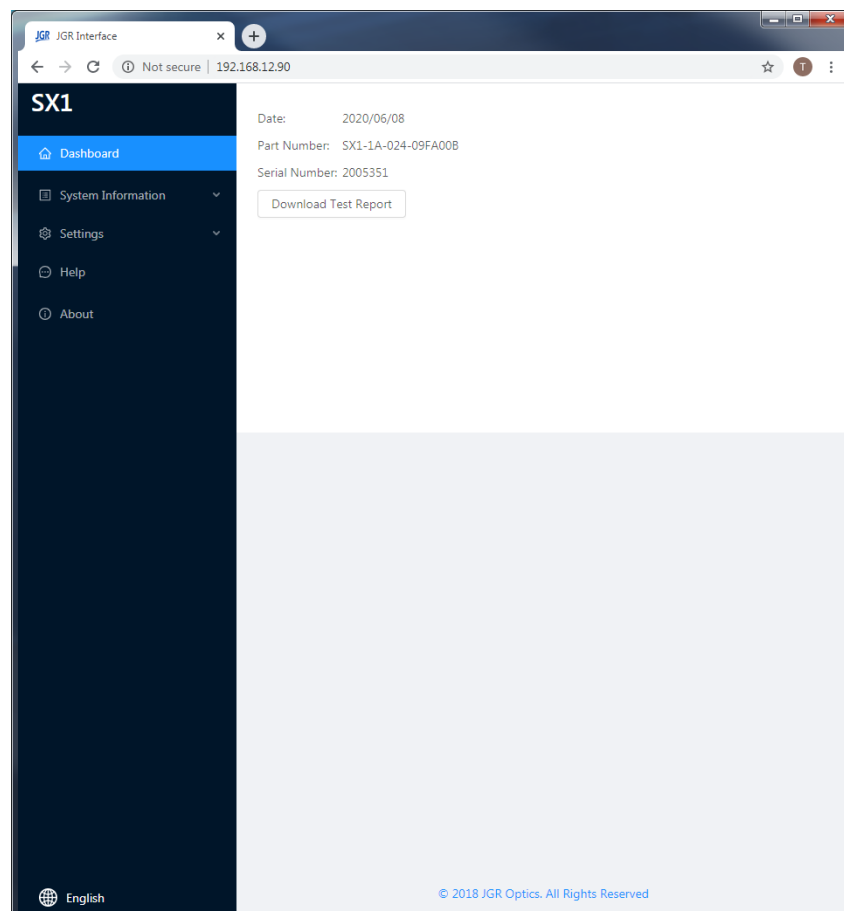


Figure 5: OSX webpage – *Dashboard* tab

Settings

Click on the *Settings* tab to expand.

Network Settings

You can view, edit or reset the network settings of the OSX from the *Settings > Network Settings* tab (Figure 6).

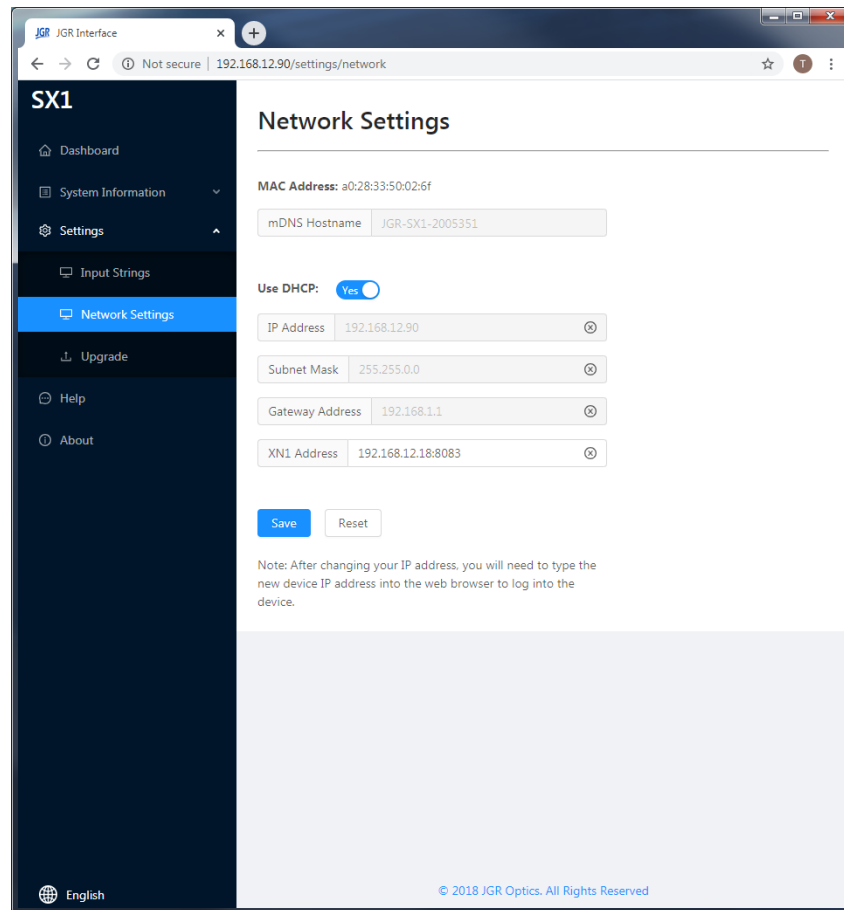


Figure 6: OSX webpage – Network Settings tab

Upgrade

Go to *Settings > Upgrade* to view the version of, upgrade or re-install the firmware of the OSX (Figure 7). Please contact support@jgroptics.com before performing a firmware upgrade for additional instructions.

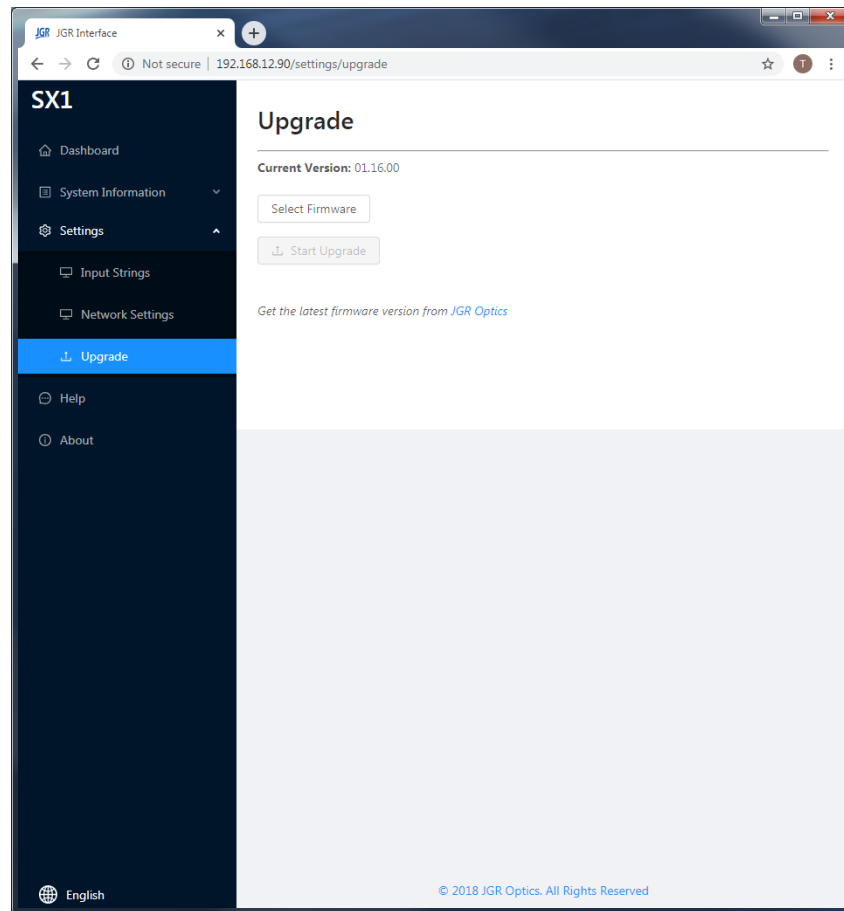


Figure 7: OSX webpage – Upgrade tab

Help

Click on the *Help* tab for the website technical support and sales contact forms as well as a link to the software downloads page.

Alternatively, email support@jgroptics.com for technical support or info@jgroptics.com for sales and other inquiries.

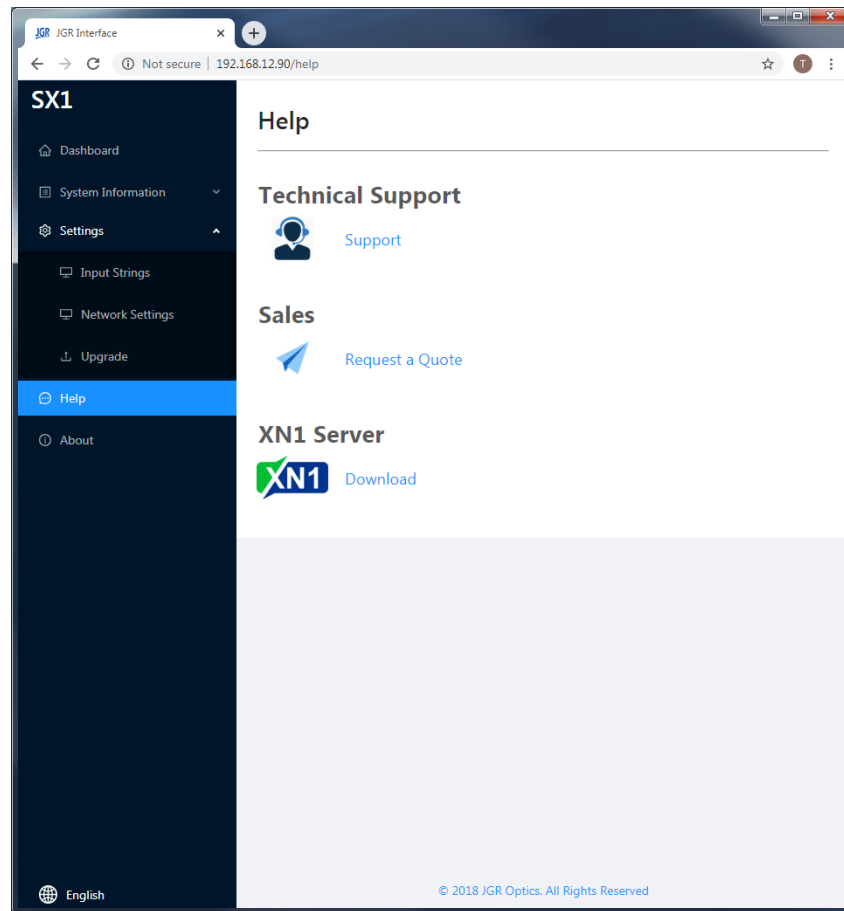


Figure 8: OSX webpage – *Help* tab

About

The *About* tab displays the unit's firmware version, model and serial number. *Advanced* mode is reserved for JGR technicians and JGR-approved service centers.

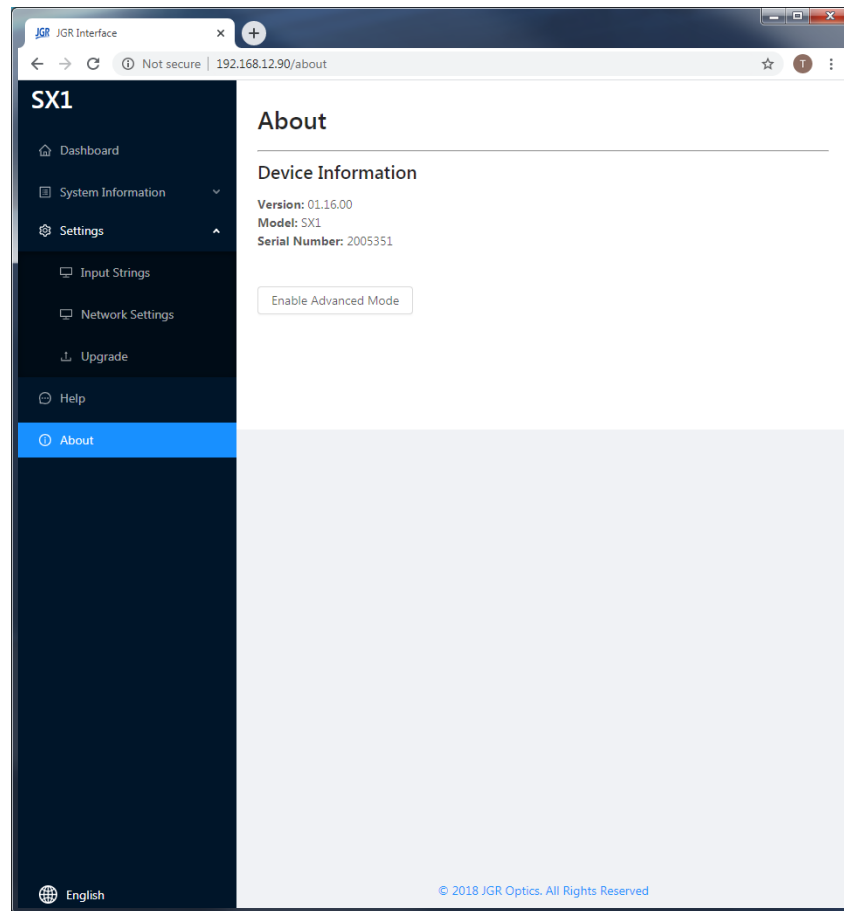


Figure 9: OSX webpage – *About* tab

7

Programming Guide

General Information

The OSX follows the *SCPI* (Standard Commands for Programmable Instruments) message-based programming standard. It conforms to the *USBTCM* (USB Test and Measurement) standard.

The TCP/IP libraries provided by most operating systems are sufficient.

Note: any *VISA* implementation is capable of controlling the OSX via TCP/IP on port 5025.

See the *Establishing Communication* section on page 10 for more information.

Step-by-step Guide

This section will provide a step-by-step programming guide in a *.NET* programming environment such as *C#* or *VB.NET*.

1. Install *VISA* drivers on the development system
2. Connect the OSX via its USB B port to the development system
3. Add a reference to *Ivi.Visa.dll* in your project:

C:\Program Files (x86)\IVI Foundation\VISA\Microsoft.NET\Framework32\v2.0.50727\VISA.NET Shared Components 5.11.0\Ivi.Visa.dll

4. Use the *IVI.Visa.GlobalResourceManager* to find all USB instruments on your system:

```
Public Overrides Function GetAllAddresses () As String ()
    Try
        Dim nameList As New List(Of String)
        nameList = GlobalResourceManager.Find("USB?*INSTR")
        Return nameList.ToArray()
    Catch ex As Exception
        Return Nothing
    End Try
End Function
```

5. Open an *IMessageBasedSession* to the desired device using an address from the *nameList* in the previous step:

```
Private visa As IMessageBasedSession

visa = GlobalResourceManager.Open(addr)
```

6. Use the *Write* method to send SCPI commands and the *Read* method to retrieve results:

```
Public Overrides Function Read(ByVal readableOnly As Boolean) As String
    Dim response As String = String.Empty
    response = visa.RawIO.ReadString()
    If response = String.Empty Then
        Throw New Exception("Read from device failed")
    End If
    Return response
End Function

Public Overrides Sub Write(ByVal strCommand As String)
    visa.RawIO.Write(strCommand)
End Sub
```

Write commands require termination with the linefeed character `\n`.

Note

The OSX runs SCPI commands asynchronously. To check if an operation is complete, it is required to poll the status bit via the query “*STAT:OPER:COND?*”. If the return value is 0, the OSX has completed its operation.

```
Try
    Write("Close" + cCommand.SChannel.ToString + vbCrLf)

    While (iResponse <> 0)
        strSwitchResponse = Query(":STAT:OPER:COND?" + vbCrLf)

        If Not Integer.TryParse(strSwitchResponse, iResponse) Then
            iResponse = -1
        End if

        System.Threading.Thread.Sleep(10)

        iCycles += 1

        If iCycles > 300 Then
            Throw New System.TimeoutException("Switch did not complete its
operation")
        End If
    End While
Catch ex As Exception
    MsgBox(ex.Message)
End Try
```

Commands Lists

See Table 4 and Table 5 for SCPI required and OSX commands respectively.

Table 4: Standard SCPI required commands list

*CLS
*ESE #
*ESE?
*ESR?
*IDN?
*OPC
*OPC?
*OPT?
*RCL "filename"
*RST
*SAV "filename"
*SRE #
*SRE?
*STB?
*TST?
*WAI
:STATus:OPERation:CONDition?
:STATus:OPERation:ENABle <byte>
:STATus:OPERation:ENABle?
:STATus:OPERation[:EVENT]?
:STATus:QUEStionable:CONDition?
:STATus:QUEStionable:ENABle <byte>
:STATus:QUEStionable:ENABle?
:STATus:QUEStionable[:EVENT]?
:STATus:PRESet
:SYSTem:ERRor[:NEXT]?
:SYSTem:VERSion?
:SYSTem:COMMunicate:LAN:ADDRes <ip/dhcp>
:SYSTem:COMMunicate:LAN:ADDRes?
:SYSTem:COMMunicate:LAN:GATEway <gateway>
:SYSTem:COMMunicate:LAN:GATEway?
:SYSTem:COMMunicate:LAN:MASK <netmask>
:SYSTem:COMMunicate:LAN:MASK?
:SYSTem:COMMunicate:LAN:HOSTname <hostname>
:SYSTem:COMMunicate:LAN:HOSTname?
:SYSTem:COMMunicate:LAN:MAC?

Table 5: OSX commands list

Command	Description
CLOSe	Move switch to the next channel.
CLOSe #	Move switch to channel #.

CLOSe?	Return current switch channel.
CFG:SWT:END?	Return switch channel count.
MODUle:CATalog?	Returns a comma separated list of all modules present in the unit. Ex. "SX 1Ax24","SX 2Bx12" List is always returned in order. I.e first item is module 1, last item is module n.
MODUle#:INFO?	Return details about module n. This info will be more detailed, equivalent to the data shown in the module selection grid. First module is Module 0.
MODUle:NUMber?	Return total number of modules present.
MODUle:SElect	Select next module.
MODUle:SElect #	Select module #.
MODUle:SElect?	Return current selected instrument.
ROUTe:CHANnel:ALL #	Change channel to # on all modules.
ROUTe#:CHANnel #	Change channel to # on module #.
ROUTe#:CLOSe #	Change channel to # on module #.
ROUTe#:CHANnel?	Return current channel of module #.
ROUTe#:COMMon #	Change common (input) on module #.
ROUTe#:COMMon?	Return currently selected common.
ROUTe#:HOMe	Reset the switch motors by finding the "home" position again. Use if switch channels ever desynchronize.
LCL #	Local = 1 or remote=0.
LCL?	Return interaction mode.
TEST:NOTIFY# "<string>"	Push a notification to the OSX LCD. # determines icon

8

Maintenance

Cleaning the Unit

1. Unplug the unit from the line power.
2. Clean the enclosure with a damp cloth.
3. Do not plug the unit back in until it is completely dry.

Cleaning the Output Connectors

1. Use a bulkhead inspection probe to inspect the output connectors.
2. If needed, use the appropriate in-bulkhead cleaning device.

9

Storage and Shipping

Damage can occur from improper handling. Make sure to maintain the unit within the specified temperature range during storage or shipping. Please follow the recommendations below to minimize the possibility of damage:

- If possible, pack the unit in its original packing material when shipping.
- Avoid high humidity or large temperature fluctuations that could generate condensation within the unit.
- Avoid unnecessary shocks and vibrations.

Returning Instruments to JGR Optics

As indicated above, please ship the returned material in the original shipping box and packing material. If these are not available, follow the guidelines below:

1. Contact JGR Optics to obtain an RMA number.
2. Cover the front panel with foam to prevent damage.
3. Wrap the unit in anti-static packaging. Use anti-static connector covers.
4. Pack the unit in a strong enough shipping box considering the unit's weight.
5. Use enough shock-absorbing material (10 to 15 cm) to cushion the unit and prevent it from moving inside the box. Pink poly anti-static foam is recommended.
6. Seal the shipping box securely.
7. Clearly mark FRAGILE on at least 3 of the 4 sides of the box.
8. Always provide the model and serial number of the unit and, if applicable, the RMA number on any accompanying documentation. If possible, indicate the RMA number on the box itself to facilitate identification.

Contact Information

JGR Optics Inc.
160 Michael Cowpland Drive
Ottawa, Ontario, Canada
K2M 1P6

Phone: +1-613-599-1000
Fax: +1-613-599-1099
Email: info@jgroptics.com
Website : www.jgroptics.com

10

Specifications

Parameter	Specification	
	Single-mode	Multimode
Wavelength Range (nm)	1250 to 1670	840 to 1350
Insertion Loss (dB) ¹	< 0.7	
Backreflection (dB) ¹	< -60	< -40
PDL (dB)	< 0.05	N/A
Repeatability (dB) ²	± 0.005	
Crosstalk (dB)	< -80	
Maximum Input Power (mW)	300	
Switching Time (ms)	300	
Control	USB or Ethernet	
Input Voltage	100 to 240 V AC, 50 to 60 Hz	
Power Consumption (VA)	60 maximum	
Switch Life (cycles)	10 ⁸	

Notes:

¹ excluding connectors

² sequential switching, add ± 0.02 dB for random

Table 6: OSX mechanical and environmental specifications sheet

Parameter	Specification	
	2 – 48 Channels	> 48 Channels
Unit Dimensions W x H x D (cm)	23.5 x 12 x 32.5 (2U half rack)	48.5 x 13.5 x 49.5 (3U full rack)
Shipping Box Dimensions W x H x D (cm)	36.5 x 39 x 53	53 x 32 x 57
Unit Weight (kg)	8	14
Total Shipment Weight (kg)	9	15
Operating Temperature (°C)	0 to 55	
Storage Temperature (°C)	-40 to 70	
Humidity (Non-condensing)	Max 95% RH from 0 to 40°C	



In the event of any trouble with this product, turn the unit off in accordance with the procedures to shut off the power described in this operation manual, disconnect the power source cord, make sure the product name and serial number described on the name plate of the product, and then contact our dealer at your place or directly contact us at Santec Photonics Laboratories. Our telephone number and facsimile number are shown below. However, we are not responsible for any trouble arising from your own repair or modification on this product.

SANTEC CORPORATION

5823 Ohkusa-Nenjyozaka, Komaki, Aichi, 485-0802, JAPAN
Tel. +81-568-79-1959
Fax +81-568-79-1718

SANTEC U.S.A. CORPORATION

Continental Plaza II, 433 Hackensack Ave., Hackensack, NJ, 07601, U.S.A.
Tel. +1-201-488-5505 (santec-1)
Fax +1-201-488-7702

SANTEC EUROPE LIMITED

Magdalen Centre, Robert Robinson Ave., The Oxford Science Park, Oxford OX4 4GA, U.K.
Tel. +44-1865-784960
Fax +44-1865-784961

SANTEC (SHANGHAI) CORPORATION, LIMITED

Shanghai Industrial Investment Building 18C, 18 Cao Xi Road (N), Shanghai 200030, P.R. China
Tel. +86-21-64279089
Fax +86-21-64279087

SANTEC CALIFORNIA CORPORATION

4750 Calle Quetzal, Camarillo, CA, 93012, U.S.A.
Tel. +1-805-987-1700

SANTEC CANADA CORPORATION

160 Michael Cowpland Drive, Ottawa, Canada, K2M 1P6
Tel. +1-613-599-1000

www.santec.com email:info@santec.com